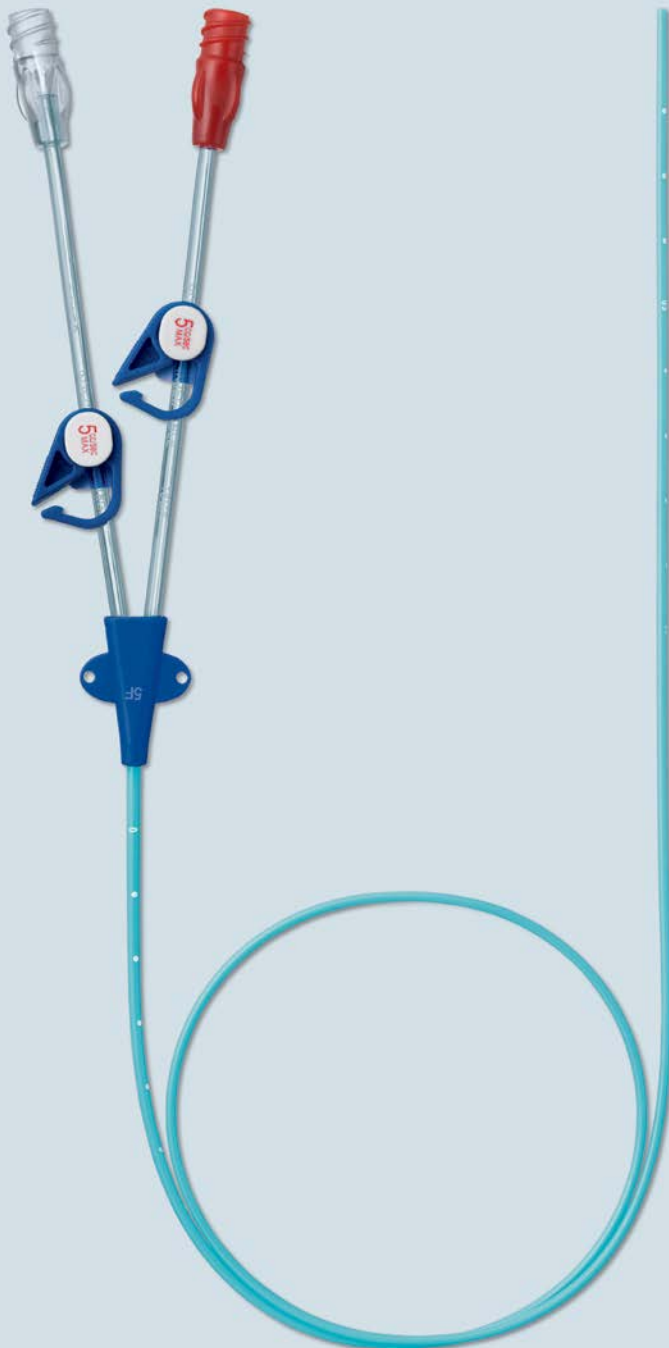

**Peripherally Inserted
Central venous Catheter
pfm CT PICC**

The pfm CT PICC is a peripherally inserted central venous catheter. It provides patient reliability, advanced technological features and comfort.

www.pfmmedical.com



The **pfm CT PICC** is a peripherally inserted central venous catheter. It provides patient reliability, advanced technological features and comfort.

Benefits

Easy identification

The blue catheter material reduces the risk of accidental misconnections to purple coloured enteral feeding devices.

Suitable for high-pressure applications

The pfm CT PICC is power injectable with 5ml / sec up to a pressure of 300 psi (21 bar).

Catheter stability

The reverse taper resists kinking and enhances catheter strength.

Various configurations

The pfm CT PICCs are available in Single-, Double- and Triple-Lumen configurations to meet the patients needs.

Greater flow rates

The increased lumen diameter has been designed to provide greater flow rates at lower pressure and easier guidewire insertion. It has been developed to minimize risk of catheter occlusion.

Details

Easy-Grip Luer-hubs

The Luer-hubs are crack resistant and combine comfort and strength. They are coloured coded for ease of identification.

Identification tags on extensions

The printed tags clearly denote maximum flow rates for contrast injection.

Technical Data

Set components:

- ▶ Blue CT-rated catheter
- ▶ Stylet 75 cm
- ▶ Tuohy-Borst Adapter
- ▶ Echogenic introducer needle 21 G / 7 cm
- ▶ Nitinol guidewire 80 cm
- ▶ Peelable sheath introducer
- ▶ Securement device
- ▶ Safety scalpel
- ▶ Syringe 10 ml
- ▶ Tape measure
- ▶ Dead end cap

Ordering Information

REF	Size	Lumen	Catheter length	Peelable sheath introducer	PU
P09070001	4.0 F (1.3 mm)	1	50 cm	4.5 F (1.5 mm)	5
P09070002	5.0 F (1.7 mm)	1	55 cm	5.5 F (1.8 mm)	5
P09070003	5.0 F (1.7 mm)	2	55 cm	5.5 F (1.8 mm)	5
P09070004	6.0 F (2.0 mm)	2	60 cm	6.5 F (2.2 mm)	5
P09070005	6.0 F (2.0 mm)	3	60 cm	6.5 F (2.2 mm)	5

Contact

Should you have any questions our Customer Solutions Team is happy to advise you.

service@pfmmedical.com
 T +49 (0)2236 9641-220
 F +49 (0)2236 9641-51

www.pfmmedical.com - Your information source for the products of pfm medical ag



pfm medical ag
 Wankelstraße 60
 50996 Köln, Germany

Certified by
 DIN EN ISO 13485

Manufacturer:
 pfm medical usa Inc.
 1815 Aston Ave, Suite 106
 Carlsbad, CA 92008, USA